

The most flexible tube rack designed for use with a multitude of tube diameters and heights

Holds virtually any brand tube 10-30mm of varying heights

Multiple Configurations

- Interchangeable tube rack inserts hold multiple tube sizes firmly in place to assist technicians constantly alternating between sample vessels
- Inserts can be pivoted in either an upper or lower position to securely accommodate the size and height of tube being used

Multi-Functional Design

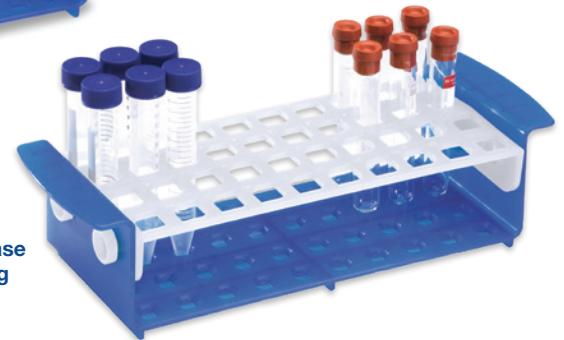
- Open design allows for better filling visibility and tube graduations
- Rack wells have a unique teardrop design that enables better retention of tubes over traditional designs
- Rotate tube rack insert to an angle position to facilitate more comfortable pipetting and allow better visibility



Lower position is ideal height when working with shorter microcentrifuge tubes

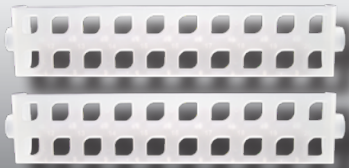
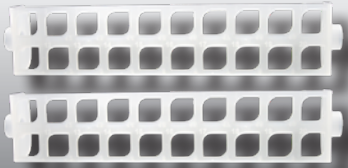



Ergonomic handles for ease of transporting the rack





Alphanumeric grid for quick sample identification and retrieval

Angled position for more comfortable pipetting

Included: 1 Base + 6 Tube Rack Inserts			
	Two 10-13 mm Tube Racks Fits 5 mL, 7.5 mL, 10 mL Tubes	Two 13-17 mm Tube Racks Fits 5 mL, 15 mL Tubes	Two 25-30 mm Tube Racks Fits 25 mL, 50 mL Tubes

OneRack® Multi

US Patent No. 9776188

Item No.	Color	L x W x H in	L x W x H cm	# Inserts Included	Well Dia/Insert	Wells/Insert	UOM
120817	Blue/Natural 	12.2 x 4.2 x 3	30.9 x 10.7 x 7.7	6 (2 ea/ Diameter)	13 mm	20	1 ea
					17 mm	20	
					30 mm	6	
120818	Purple/Natural 	12.2 x 4.2 x 3	30.9 x 10.7 x 7.7	6 (2 ea/ Diameter)	13 mm	20	1 ea
					17 mm	20	
					30 mm	6	